



apex locators - micromotors

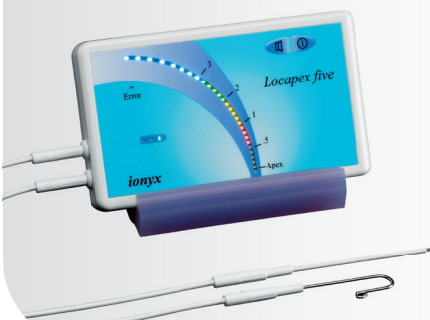



Endy 6200

micromotor + apex locator



Learn about the other products in our [ionyx](#) range

	<p>Locapex Five apex locator Apex locator with instant measurement system</p>	<p>Endy NT2 micromotor with light Portable endodontic micromotor with built-in light</p> 
--	--	---

Endy 6200

micromotor + apex locator

instant measurement system

The instant measurement system provides you with information on the position of your file in the canal the moment it makes contact with the canal.

rapid motor shutdown

As soon as the programmed torque is achieved, the motor is quickly shut down by means of an electronic lock, limiting the risk of breakage of the instrument.

a speed range suited to root canal files

The Endy 6200 lets you set your file speed between 100 and 500 revolutions per minute.

slow reverse

This 1-second slow-reverse technology used by IONYX since 1996, lets you remove the instrument without sending dental debris back towards the apex.

12-memory settings

The Endy 6200 has a 12-memory settings for programming as many torques and speeds corresponding to each of the files used. Programming is easy and straightforward.

long battery life

The Endy 6200 lets you work on battery power for 1 hour with an average torque of 1.5 N.cm (40 root canal treatments on average reported).

quick charger

The device's quick charger fully recharges your battery in 2 hours. In order to preserve your battery, the charge automatically shuts down as soon as the maximum is reached.

The Endy 6200 electronic casing has a 2-years warranty.
The contra-angle and the battery have a 1-year warranty.

Product strengths

100% digital measurement

The signal emitted by your Endy 6200 is 100% digital. There is no drift or modification over time in the accuracy of your IONYX apex locator.

a built-in insulated sensor in the contra-angle

As soon as your file is inserted into the contra-angle, the measurement sensor is automatically connected. The Peek plastic contra-angle ensures total insulation from the measurement signal.

consistent torque and speed

The speed of the Endy 6200 remains consistent regardless of the force with which the file is used in the root canal. It never varies, even if the battery is less than 50% charged. The torque is always controlled in order to preserve your files.



*Your partner
for over 20 years*

ionyx

34 rue de Fleurence
F-33290 BLANQUEFORT

T. +33 5 56 35 34 27
F. +33 5 56 95 16 20

ionyx@ionyx.eu
www.ionyx.eu

ionyx
endodontic
treatment